MEDISANA®

home of wellness

DE Shiatsu-Massagegerät NM 860

GB Shiatsu Massager NM 860

FR Appareil de massage Shiatsu NM 860

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IMPORTANT INFORMATION! RETAIN FOR FUTURE USE!

Read the instruction manual carefully before using this device, especially the safety instructions, and keep the instruction manual for future use.

Should you give this device to another person, it is vital that you also pass on these instructions for use.

Explanation of symbols

This instruction manual belongs to this device. It contains important information about starting up and operation. Read the instruction manual thoroughly. Non-observance of these instructions can result in serious injury or damage to the device.

WARNING

These warning notes must be observed to prevent any injury to the user.

CAUTION

These notes must be observed to prevent any damage to the device.

NOTE

These notes give you useful additional information on the installation or operation.

Protection category II

LOT

Lot number

Manufacturer

about the power supply

- Before connecting the device to your power supply, please ensure that the supply voltage stated on the rating plate is compatible with your mains supply.
- Only insert the power supply unit into a socket when the appliance is switched off.
- Only use the power supply unit supplied for the massager.
- Keep the power supply unit, cable and appliance well away from heat, hot surfaces, moisture and liquids.
- Do not touch the power supply unit when you are standing in water and always take hold of the plug with dry hands.
- Do not take hold of a device if it falls into water. Pull out the power plug and the power supply unit from the socket immediately.
- Connect the device only in that way, so that the power plug is freely accessible.
- Always switch off the appliance after use with the buttons on the control unit and unplug the power supply unit.
- To disconnect the unit from the power supply, always pull the plug out of the socket. Do not pull on the cable!
- Do not carry, pull or twist the device on the cable or the power supply unit.
- You must not continue to use the unit if the lead or the mains adaptor are damaged. For safety reasons, these parts may only be replaced by authorized service personnel. In order to prevent hazards, always send the unit in to the service centre.
- Ensure that the cable will not cause tripping. The cable must not be kinked, trapped or twisted.

specific groups of people

- This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
- Children must be supervised to ensure that they do not play with the device.
- Do not use this unit to supplement or replace medical treatment. Chronic complaints and symptoms may be made worse.
- Do not use the massager, or only after consulting your doctor, if:
 - you are pregnant,
 - you have a cardiac pacemaker, artificial joints or electronic implants,
 - you suffer from one or more of the following illnesses and complaints: Circulation disorders, varicose veins, open wounds, bruises, broken skin, inflammation of the veins.
- Do not use the appliance near your eyes or near any other sensitive areas.
- If you experience any pain or discomfort during the massage, stop immediately and consult your doctor.
- Take care when using the heating function. **The appliance has hot surfaces.** Do not use the appliance on small children, vulnerable individuals or people who are not sensitive to heat.
- Before using the massager, consult your doctor regarding any unexplained pain, if you are undergoing medical treatment or using medical equipment.

before using the unit

- Check the power supply unit, cable and massager carefully for damage before each use. Do not attempt to use a defective unit.
- Do not use the appliance if it, or the power cord show signs of damage, if it is not in perfect working order, or if it has been dropped or become damp. In order to prevent hazards, always send the unit in to the service centre.

operating the unit

The device is not to be used for commercial use or medical purposes.

If you have health concerns, consult your doctor before using the device.

Only use the massager indoors.

Do not use the massager in damp rooms (e.g. when bathing or showering).

- The device must only be used for its intended purpose as described in the instruction manual.
- Using the unit for any other purpose invalidates the warranty.
- Only use the appliance in the proper place, as described in this instruction manual.
- Do not use the appliance while you are asleep or when you are lying in bed.
- Do not use the appliance before going to sleep. The massage has a stimulating effect.
- Never use the appliance while driving a vehicle or operating a machine.
- The maximum operating time for one session is 15 minutes. Prolonged use reduces the service life of the appliance.

- Very long periods of use can lead to overheating. Always allow the appliance to cool down before using it again.
- Never put the unit down next to an electric stove or other sources of heat.
- Never leave the unit unattended if it is connected to the power supply.
- Do not allow the massage cushion to come into contact with pointed or sharp objects.
- Do not cover up the appliance when it is switched on. Never use it under blankets or cushions. There is a risk of fire, electric shock or injury.
- Make sure, that no fingers or other parts of the body come between the rotating massage heads ❶ . Risk of injuries!
- Do not use an external time switch or a separate remote control system to operate the appliance.

maintenance and cleaning

- The unit is maintenance free. If a fault does occur, just check whether the mains power adapter is properly connected.
- You may only clean the device yourself. If the device malfunctions please do not repair it yourself as this leads to the automatic expiry of all guarantee rights. Ask your specialist dealer and have repairs carried out by the service.
- Cleaning and user maintenance shall not be made by children without supervision.
- Do not immerse the device in water or other fluids.
- If, nonetheless, liquid should penetrate the device, unplug the mains cable immediately.

Do not wash! Only the cover ❷ is washable. For details, refer to chapter 4.1 *Cleaning and maintenance*.

Do not dry clean!

Thank you very much

Thank you for your confidence in us and congratulations!
You have gained a high-quality product from **MEDISANA** with the purchase of this Shiatsu Massager **NM 860**. This appliance is intended for the Shiatsu massage of the neck and shoulder area - additionally it can be flexibly used on other parts of the body (e.g. on the back, waist, upper and lower leg). To derive the most success and to enjoy your **MEDISANA** Shiatsu Massager **NM 860** for a very long time, we recommend that you read the following information on the use and care of the appliance carefully.

**2.1
Items
supplied and
packaging**

Please check first of all that the device is complete and is not damaged in any way. In case of doubt, do not use it and contact your dealer or your service centre.

The following parts are included:

- 1 **MEDISANA** Shiatsu Massager **NM 860**
- 1 Power supply with cable
- 1 Instruction manual

The packaging can be reused or recycled. Please dispose properly of any packaging material no longer required.

If you notice any transport damage during unpacking, please contact your dealer without delay.

WARNING

Please ensure that the polythene packing is kept away from the reach of children! Risk of suffocation!

**2.2
What is
Shiatsu
massage?**

Shiatsu is a form of finger pressure massage and one of the most important massaging techniques along with manual lymphatic drainage and reflexology. This form of physical therapy was developed in Japan but stems from traditional Chinese massage. Gentle stroking and beneficial pressure application techniques are used to awaken life forces and encourage them to flow through the body.

The object of the massage is not just to ease tension, but also to optimise the flow of life force and harmonise body and soul.

**2.3
How does
the
MEDISANA
Shiatsu
Massager
NM 860
work?**

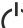
In the **MEDISANA Shiatsu Massager NM 860** you have a modern device for genuine Shiatsu style massages on different parts of the body. The **MEDISANA Shiatsu Massager NM 860** has four rotating massage heads for Shiatsu massage. The device is equipped with an additional warming function in the massage heads. The unique combination of intensive shiatsu massage with radiated heat makes the massage an effective and relaxing treatment. The settings can be comfortably adjusted by means of the easy-to-handle integrated buttons. The hand straps allow an individually adjusted massage intensity by pulling the straps accordingly.

CAUTION

Make sure that the maximum operating time of 15 minutes is not exceeded.

3.1 Neck massage

1. Insert the power supply unit plug into the connection for power supply  on the backside of the device.
2. Plug in the power supply unit .
3. Place the massager around the neck so that it has sufficient hold on the shoulders.
4. Guide your hands through the hand straps  on right and left side. By pulling on the straps you may stabilize the device in its position and adjust the massage intensity (pulling more stronger increases the massage intensity).
5. Press the ON-/OFF-button  to switch on the device. The Shiatsu massage starts on the slow setting.
6. By pressing the MASSAGE-button  you may change the massage setting to fast speed. Pressing the MASSAGE-button  again to switch the appliance off.
7. To switch on the warming function, press the WARMTH-button . The warming function only works with active massage function and generates a comfortable, smooth warmth.
8. To switch off the warming function, press the WARMTH-button  again.
9. Press the ON-/OFF-button  to switch off the device. After 15 minutes of constant massage the device will switch off automatically.

3.2 Massage of other parts of the body

The device can alternatively also be used on other parts of the body, e.g. on the back, shoulders, hips, upper and lower legs. To apply the Shiatsu massage there also, put the appliance on the decided area of your body and hold it there during the massage. You may also use the additional warming function.

4.1

Cleaning and maintenance

- Before cleaning the appliance, make sure that the appliance is switched off and the power supply unit is disconnected. Allow the unit to cool down.
- Clean the massager with a slightly damp sponge only. Never use brushes, aggressive cleaning agents, petrol, thinners or alcohol. Rub the seat dry with a clean dry cloth.
- Never immerse the appliance in water to clean it and ensure water cannot get into the appliance.
- Wait until the device is completely dry before using it again.
- Undo any kinks in the cable.
- Take care when rolling up the mains lead to avoid damaging it.
- You should ideally store the device in its original packaging and keep it in a clean, dry place.

The washable cover ② of the massage heads can be removed for washing.

You may wash it in a normal wash at max. 30 °C. Follow the notes on the label. To remove the cover ②, open the zip fastener all around it and take it off carefully.

4.2

Disposal

This product must not be disposed of together with domestic waste. All users are obliged to hand in all electrical or electronic devices, regardless of whether or not they contain toxic substances, at a municipal or commercial collection point so that they can be disposed of in an environmentally acceptable manner.

Consult your municipal authority or your dealer for information about disposal.

4.3

Technical specifications

- | | |
|----------------------|--|
| Name and model | : MEDISANA Shiatsu Massager
NM 860 |
| Voltage supply | : Input: 100 - 240V~, 50/60 Hz, 500mA
Output: 12V=1500 mA |
| Power consumption | : 18 watts |
| Automatic switch-off | : After approx. 15 minutes |
| Operating conditions | : Only use in dry rooms |
| Storage conditions | : Clean, dry and cool |
| Dimensions | : Approx. 58 x 34 x 20 cm |
| Weight | : Approx. 1,55 kg |
| Article No. | : 88942 |
| EAN Code | : 40 15588 88942 4 |

In accordance with our policy of continual product improvement, we reserve the right to make technical and optical changes without notice.

The current version of this instruction manual can be found under www.medisana.com

Warranty and repair terms Please contact your dealer or the service centre in case of a claim under the warranty. If you have to return the unit, please enclose a copy of your receipt and state what the defect is.

The following warranty terms apply:

1. The warranty period for **MEDISANA** products is three years from date of purchase. In case of a warranty claim, the date of purchase has to be proven by means of the sales receipt or invoice.
2. Defects in material or workmanship will be removed free of charge within the warranty period.
3. Repairs under warranty do not extend the warranty period either for the unit or for the replacement parts.
4. The following is excluded under the warranty:
 - a. All damage which has arisen due to improper treatment, e.g. non-observance of the user instructions.
 - b. All damage which is due to repairs or tampering by the customer or unauthorised third parties.
 - c. Damage which has arisen during transport from the manufacturer to the consumer or during transport to the service centre.
 - d. Accessories which are subject to normal wear and tear.
5. Liability for direct or indirect consequential losses caused by the unit are excluded even if the damage to the unit is accepted as a warranty claim.

MEDISANA AG

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41468 NEUSS

Germany

eMail: info@medisana.de

Internet: www.medisana.com

The service centre address is shown on the attached leaflet.

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